

NEW!



IntelliPAP AutoBilevel System
(Model DV57D-HH)

IntelliPAP® AutoBilevel System

Comfortable therapy at higher pressures

Introducing the latest addition to the IntelliPAP line — the IntelliPAP AutoBilevel system. This system **incorporates the comfort of AutoAdjust® and Bilevel technology** to deliver comfortable therapy particularly at higher prescription pressures. The IntelliPAP AutoBilevel **automatically adjusts to the user's pressure requirements** on both inhalation and exhalation, using the unique DeVilbiss AutoAdjust algorithm featuring Event Set Measurement.


Performance Features

- **Includes SmartCode® Therapy Monitoring** to remotely track adherence data
- **Optional SmartLink® Therapy Management** system available for remotely monitoring usage and effectiveness of therapy
- **Patented Flow Rounding™ Technology**
- Fixed Pressure Support – adjustable from 0 to 12 cmH₂O
- Respiratory Event Detection
- Automatic pressure adjustment
- Adjustable definition of apnea and hypopnea
- **Among the quietest AutoBilevel units on the market**
- Easy-to-read LCD display
- Pressure delay settings
- Snoring and mask condition detection
- Optional integrated heated humidification system
- Automatic Altitude Adjustment up to 8,500 feet
- **Industry-leading 3-year warranty**

Design Features

- Whisper-quiet operation – 26 dBA!*
- Exceptionally small footprint fits easily into the home environment and carry-on luggage
- Includes discreet carry bag that can be easily expanded to hold additional accessories
- Optional integrated heated humidification system easily docks into the bottom of the unit without requiring a separate hose or power cord
- Humidifier reservoir is easy to remove and replace, and includes built-in overflow protection
- Water ingress protection prevents humidifier water from damaging the CPAP
- Accommodates various tubing lengths (6 and 10 feet)
- Large, slip-resistant pad on bottom of unit to minimize movement on night table during the night
- Eight-foot power cord length provides greater positioning flexibility – no bulky power brick
- Easy and intuitive set up – three, large buttons control nearly all functions for simplified patient operation
- Large, backlit LCD and easy-to-follow patient menu
- Auto on/off** automatically starts the flow of air after breathing into the mask and stops the flow of air when mask is removed
- Delay pressure eases the machine to the prescribed pressure – allowing the user to fall asleep at a lower pressure (0-45 minutes in 5-minute increments and .5 increments for pressure)
- Reminders alert user to replace or check parts
- Extra air inlet and fine particle filters included

Selectable IPAP and EPAP Flow Rounding*** (0-5)

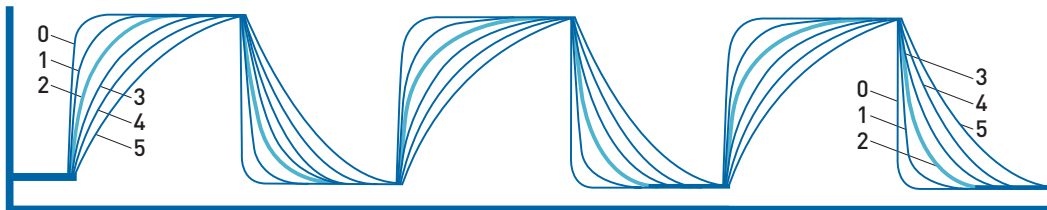

Technical Specifications

Product Dimensions	6.4" H x 6.5" W x 8.4" D	Automatic Leak Compensation	Yes
Product Weight	4.45 lbs with humidifier	Auto On/Off	Yes (can be enabled or disabled)
Pressure Range	3-25 cm H ₂ O (.5 increments)	Motor Hour Meter	Yes
Universal Power Supply	Yes	Therapy Compliance Meter	Yes
AC Cord Length	8 feet	Warranty	3-year warranty on IntelliPAP AutoBilevel 3-year warranty on humidification heater
DC Operable	Yes	DeVilbiss SmartLink Compatible	Yes
Automatic Altitude Adjustment	Yes	SmartCode Feature	Yes
ABOB Analog Output for PSG	Yes		

* Tested per ISO 17510

** U.S. Patent 5,503,146

*** U.S. Patent 5,865,173



DeVilbiss Healthcare LLC | 100 DeVilbiss Drive | Somerset, PA 15501 USA
800.338.1988 | 814.443.4881 | www.DeVilbissHealthcare.com