

Nonin[®] Xpod[®] Oximeter

A Powerful, Clinical Option for Superior, Patient Therapy Management

The Nonin Xpod Pulse Oximeter can now be paired as an optional accessory with the DeVilbiss SmartLink[®] Therapy Management System to provide up to nine nights* of patient oximetry data.

While SmartLink monitors and records up to 72 hours of continuous SpO₂ and pulse data from the Xpod Oximeter, the IntelliPAP[®] Standard, Standard Plus, AutoAdjust[®], or Bilevel model displays the information, which can also be downloaded to SmartLink PC software and viewed in report format.

Included in this accessory package is the Nonin soft sensor finger probe that comfortably and securely fits onto a user's finger; it stays in place for overnight use. SmartLink records only when the oximeter receives a strong signal — no more interference from ambient light.

Performance Features

- SpO₂: range 0%-100%; accuracy — +/-1% at range 70%-100%
- Pulse rate: range 18-300 bpm; accuracy — +/-3 bpm
- Sample rate: 0.25 Hz (SpO₂ and pulse rate)




* 9 nights x 8 hours/night = 72 hours total recording time



- PureSAT® technology provides exceptionally accurate pulse oximetry
- A full range of sensors can be used with the Xpod — reusable, flex, and soft sensor oximeter
- The Xpod Oximeter delivers superior performance and can be used in motion and low perfusion
- The Xpod Oximeter uses a low power draw, only 60 mW

Design Features

- The soft sensor finger probe is comfortable and stays in place, making it perfect for overnight use
- The Xpod Oximeter and soft sensor are designed to integrate with SmartLink, using only the Xpod cable
- Form-fitted design decreases ambient light interference
- Durable and easy to clean

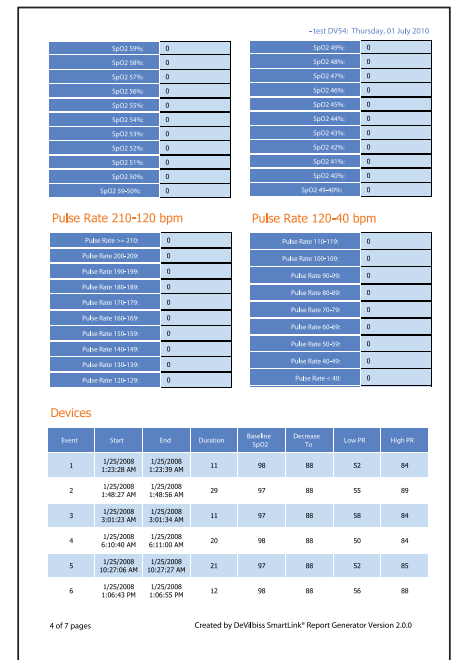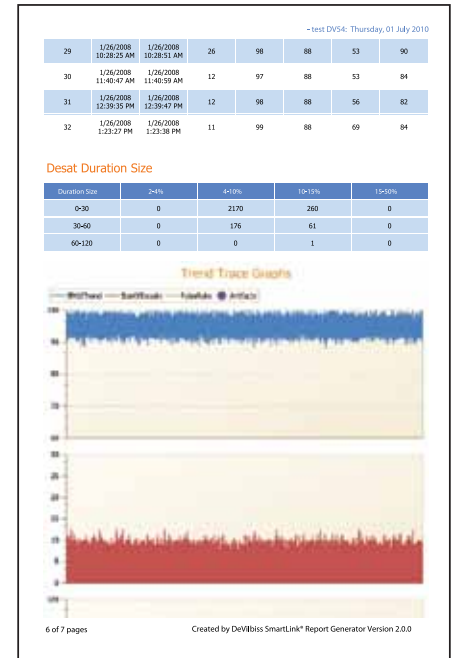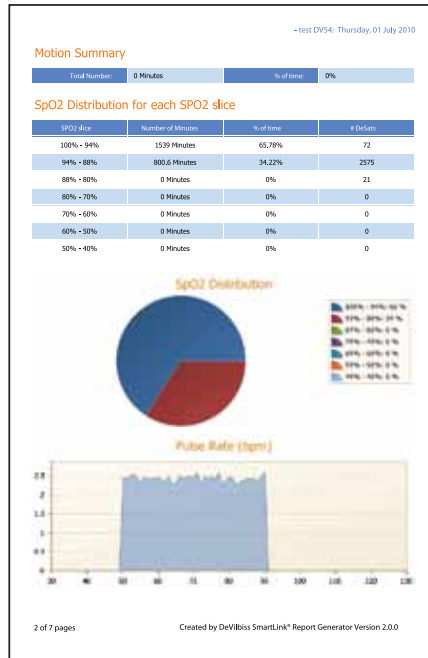

Accessories and Replacement Parts

Part Number	Product Description	Part Number	Product Description
DV5M-697	Oximeter Kit	9054D-695	Flex Oximeter Probe, Reusable
8000AA	Oximeter Finger Clip, Reusable	8000SM	Oximeter Finger Soft Sensor, Medium
8000Q	Oximeter Ear Clip, Reusable	8000SL	Oximeter Finger Soft Sensor, Large
9054D-693	Oximeter Patient Cable	8000JFW	FlexiWrap® Disposable Adhesives (25 pack)



DeVilbiss Healthcare LLC | 100 DeVilbiss Drive | Somerset, PA 15501 USA
 800.338.1988 | 814.443.4881 | www.DeVilbissHealthcare.com